

Diagnopal Dynabioassay Products LTD

Diagnopal Interference Binder (Diluent-Reaction Buffer)

DPIBA (100ml), DPIBA (1000ml)

Introduction:


Diagnopal Interference Binder/Diluent is an assay buffer for universal Heterophilic Antibody neutralizing reaction matrix that eliminates a great range of noises (HAMA, RF, HF) in Immunoassay tests. It Blocks Non-Specific Binding of sample analyte and antibodies as well and makes significant increase in signal to noise ratio in immunoassay tests. It is recommended to be used in ELISA, CLIA, Lateral flow, RIA and Western blotting. Just easily dilute your samples in 1/2 to 1/100 with Diagnopal Interference Binder/Assay Buffer and perform your test according to your own standard protocol or test method!

Specification

Stability and storage	The product should be kept at 2-8 °C. In refrigerated condition it is stable for 3 years from date of manufacture.
Specification	Protein: BSA/others
	Product Buffer: PBS
	pH: 6.4-7.2
	Preservative: Methyl Isothiazolone
	Acceptable Diluents: Do not dilute
	Dry weight(gr/mL): 0.048-0.064
Shipment Temperature: RT	
Notes	Keep product away from freezing, Direct sun light and heat. 3 months of storage of un-opened products in Dark and RT will not affect the product quality.
Available forms	Ready to use

Diagnopal Dynabioassay Products LTD

Example of Diagnopal Interference Binder/Assay Buffer performance in real test

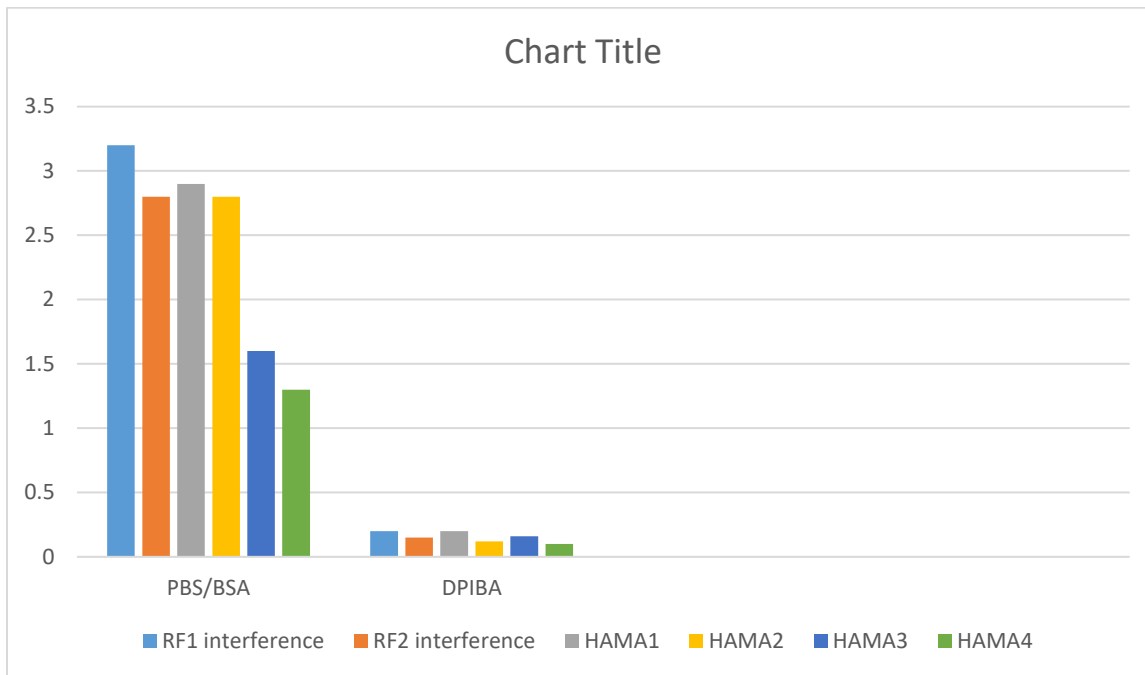

Recommendation for use:

Consumer can directly dilute sample/Serum (1/2-1/100) in Diagnopal Interference Binder/Assay Buffer for ELISA, CLIA, RIA tests.

This product is for research and/or further manufacture use only and is not for human therapeutic or consume purpose. Sales without any seller's warranty or representation. No claims for replacement or refund of purchase shall be allowed. All claims must be made within 30 following date of delivery.